

# Premium Patella Saw Guide


Providing excellent stability and control during patella resection.

Featuring an ergonomic pistol grip the Enztec Premium Patella Saw Guide allows medial-lateral clamping, protecting the patella tendon and giving you confidence when making your cuts.

The Premium Patella Saw Guide is available in three variants; 9mm and 10mm base resection heights (compatible with removable shims) and an 8-10mm fixed stylus version. The tooth profile provides excellent patella engagement and the dual lock buttons prevent accidental release.

Loved by surgeons worldwide, you can rely on the Enztec Premium Patella Saw Guide to keep you in control.



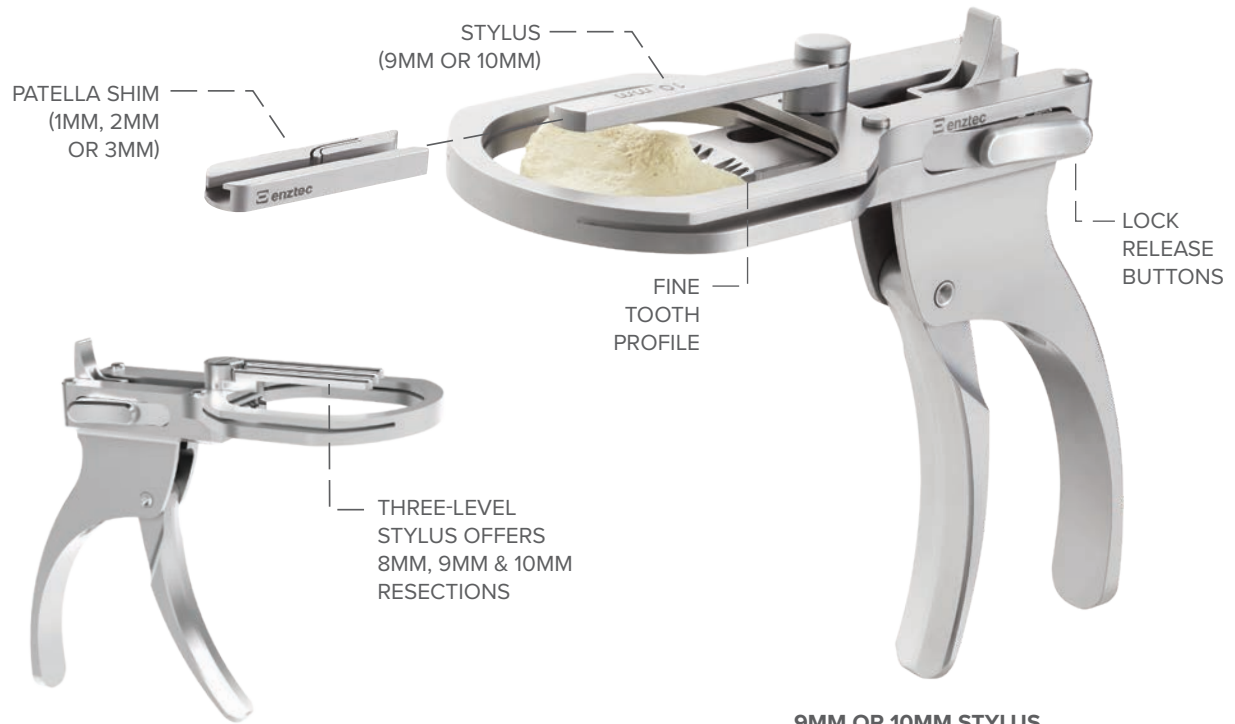

**8-10MM FIXED STYLUS**

**9MM OR 10MM STYLUS**

## Features

- Three variants are available (9mm, 10mm and 8-10mm fixed stylus)
- Three removable shims (-1mm, -2mm or -3mm) for resection height adjustment
- For use with a 1.27mm or 1.37mm saw blade
- Wide excursion, dual sided saw slot
- Rapid clamping and locking mechanism

## Specifications

- Materials: Stainless Steel
- Weight: 630g / 1.4lbs
- Length: 175mm / 6.8in
- Add your logo: available