




Modular Offset Cup Impactor X/XR

Modular design for workflow flexibility and visualisation.

The Enztec Modular Offset Cup Impactor combines trusted design with added versatility. Configured to connect to powered devices or Enztec's Modular Anvil, it allows both powered and manual impaction, giving you flexibility in any workflow.

Its detachable trinket allows quick in-situ visualisation or intraoperative imaging to verify cup position before final impaction. The generous offset supports a variety of patient anatomies, making this instrument suitable for any approach.

Available in two variants, Enztec Modular OCI X for standard threaded cups and the Enztec Modular OCI XR for cups with an anti-rotation feature.

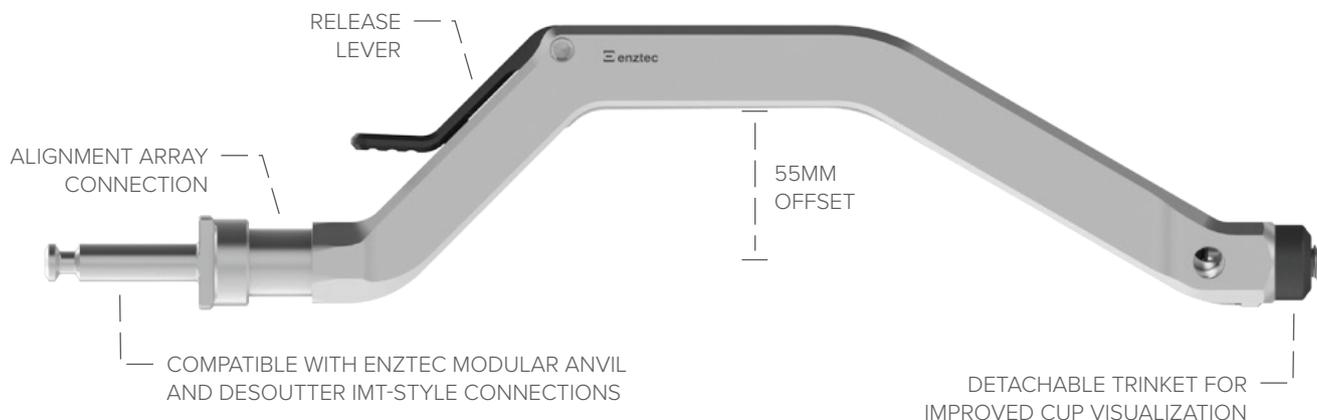

COMPATIBLE IMPLANTS


MODULAR OCI X

- S+N R3 & REDAPT
- DePuy Pinnacle & Emphasys
- Stryker Trident
- Signature Logical G
- Microport Procotyl P
- Corin Trinity
- + more

MODULAR OCI XR

- Zimmer Biomet G7
- Globus Onvoy
- TJO Klassic
- + more

TRINKET TYPES

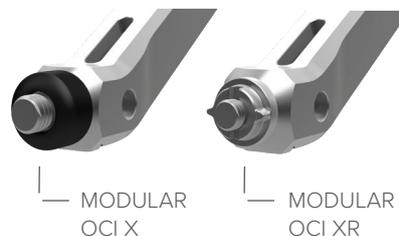

Features

- Compatible with DeSoutter IMT Connection and Enztec Modular Anvil for versatility between powered and manual impaction
- 55mm offset for a wide range of patient anatomies
- Detachable trinket for quick visualisation or intraoperative imaging
- Easy release lever for straightforward handling
- Lock mechanism to set cup orientation prior to impaction
- Alignment array connection for accurate placement

Specifications

- Materials: Stainless Steel, Silicone, Propylux
- Weight: 722g / 1.6lbs
- Length: 325mm / 12.8in
- Offset: 55mm / 2.2in
- Add your logo: available