




Partial Knee Impactor with Slap Hammer

Providing femoral trial extraction as well as impaction of final components.

The Enztec Partial Knee Impactor builds on the proven design of the Enztec Total Knee Femoral Impactor—one of our most recognised instruments, with over 15 years of clinical history.

Designed specifically for PKA/UKA, it features a polymer impaction head with three points of contact giving surgeons greater control during guided placement and final impaction.

This variant includes a slap hammer handle, allowing for easy femoral trial removal while maintaining full functionality for final impaction and secure implant seating.

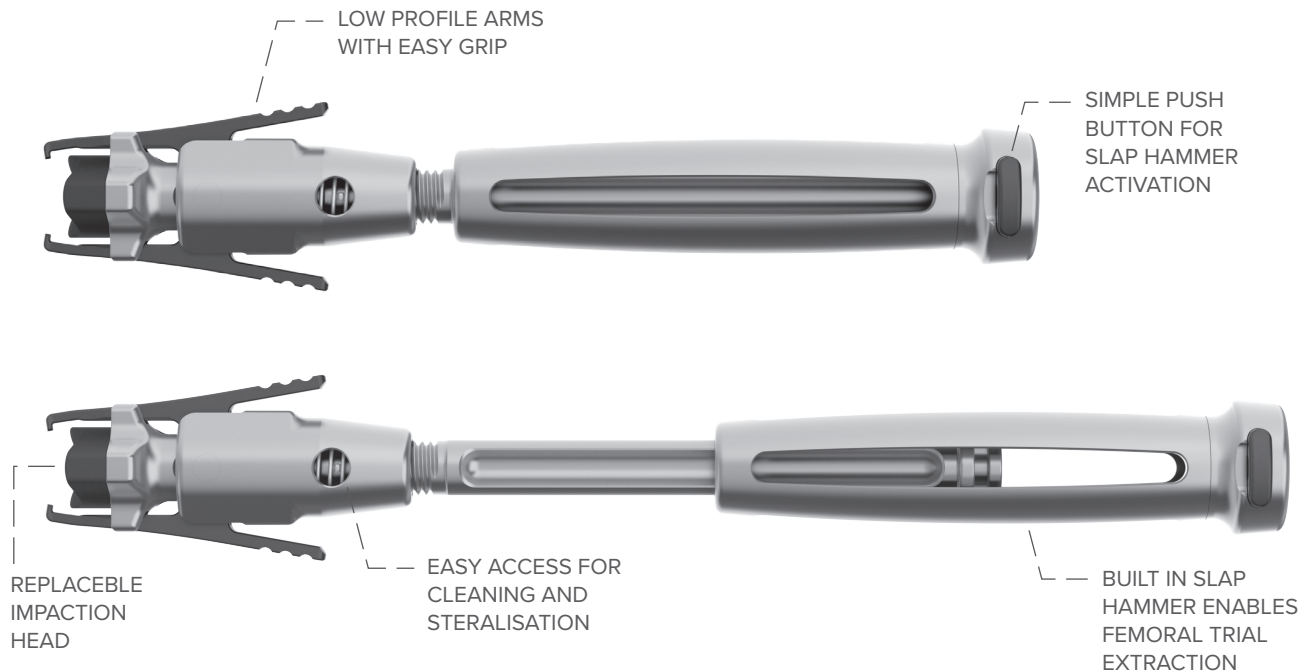

Features

- Intuitive to use
- Integral slap hammer handle enables seamless femoral trial extraction
- Configurable to suit your femoral components
- Open design ensures it is easy to clean
- Replaceable impaction head.

Specifications

- Materials: Stainless steel, Propylux
- Weight: 428g / 0.9lb
- Length: 215mm / 8.4in
- Add your logo: available