



Classic Offset Cup Impactor X/XR


A new level of in-situ visualization.

With its detachable trinket the Enztec Classic Offset Cup Impactor can be quickly and easily removed for visual inspection in-situ or intraoperative imaging to verify cup position before final impaction.

The generous offset supports a variety of patient anatomies meaning this instrument can be used in any approach. It is already compatible with a wide range of acetabular cups or will be adapted to fit your cup.

The Offset Cup Impactor has two variants; the OCI X for standard threaded cups and the OCI XR for threaded cups containing an anti-rotation feature. Giving you flexibility and confidence in any situation.



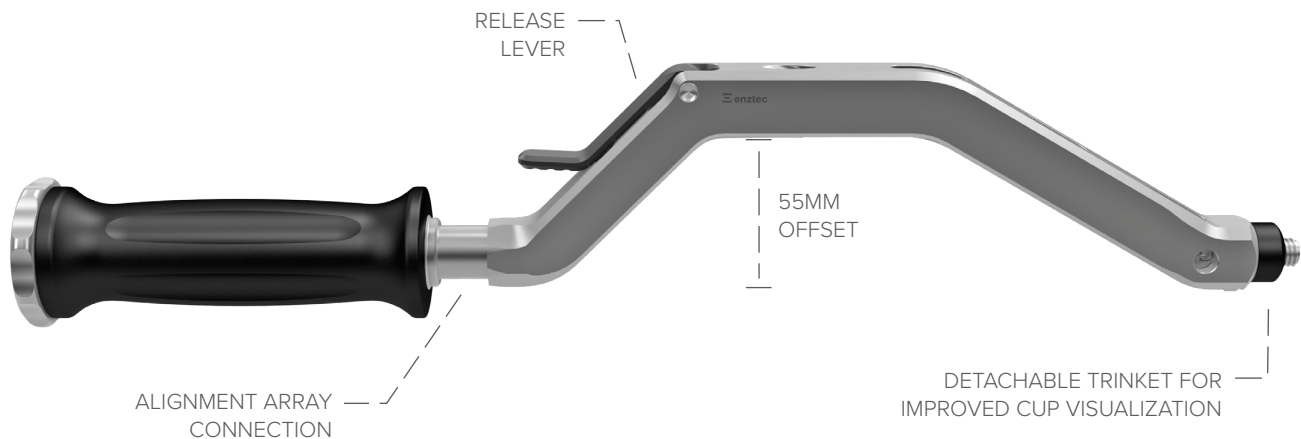

COMPATIBLE IMPLANTS


OCI X

- S+N R3 & REDAPT
- DePuy Pinnacle & Emphasys
- Stryker Trident
- Signature Logical G
- Microport Procotyl P
- Corin Trinity
- + more

OCI XR

- Zimmer Biomet G7
- Globus Onvoy
- TJO Klassic
- + more

TRINKET TYPES

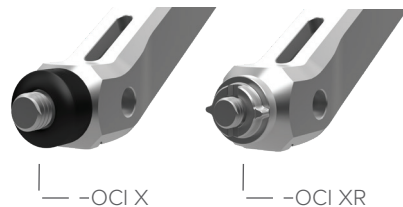

Features

- Allows in-situ visualization and intraoperative imaging of the cup before final impaction
- Lock mechanism to set the cup orientation prior to impaction
- Release lever clear of the incision
- Compatible with Enztec's Alignment Arrays and NAVBIT SPRINT™

Specifications

- Materials: Stainless Steel, Propylux, Silicone
- Weight: 1300g / 2.9lbs
- Length: 420mm / 16.5in
- Offset: 55mm / 2.2in
- Add your logo: available