


Premium Patella Saw Guide

Premium Patella Saw Guide

The Premium Patella Saw Guide provides a fast and ergonomic method of gripping the patella for resection.

Features:

- Two variants are available offering either a 9mm or 10mm standard resection
- Three shims are available (1mm, 2mm, 3mm), which allow for bone preservation, reducing the height of the resection by up to 3mm
- Fine tooth profile with mixture of long/short teeth provides positive engagement with patella
- Pistol grip for ergonomic handling
- Dual release lock buttons prevent accidental release
- Finger grip located on jaw allows for forceful release, should teeth stick in patella
- Simple proportional movement of clamping jaw relative to handle.

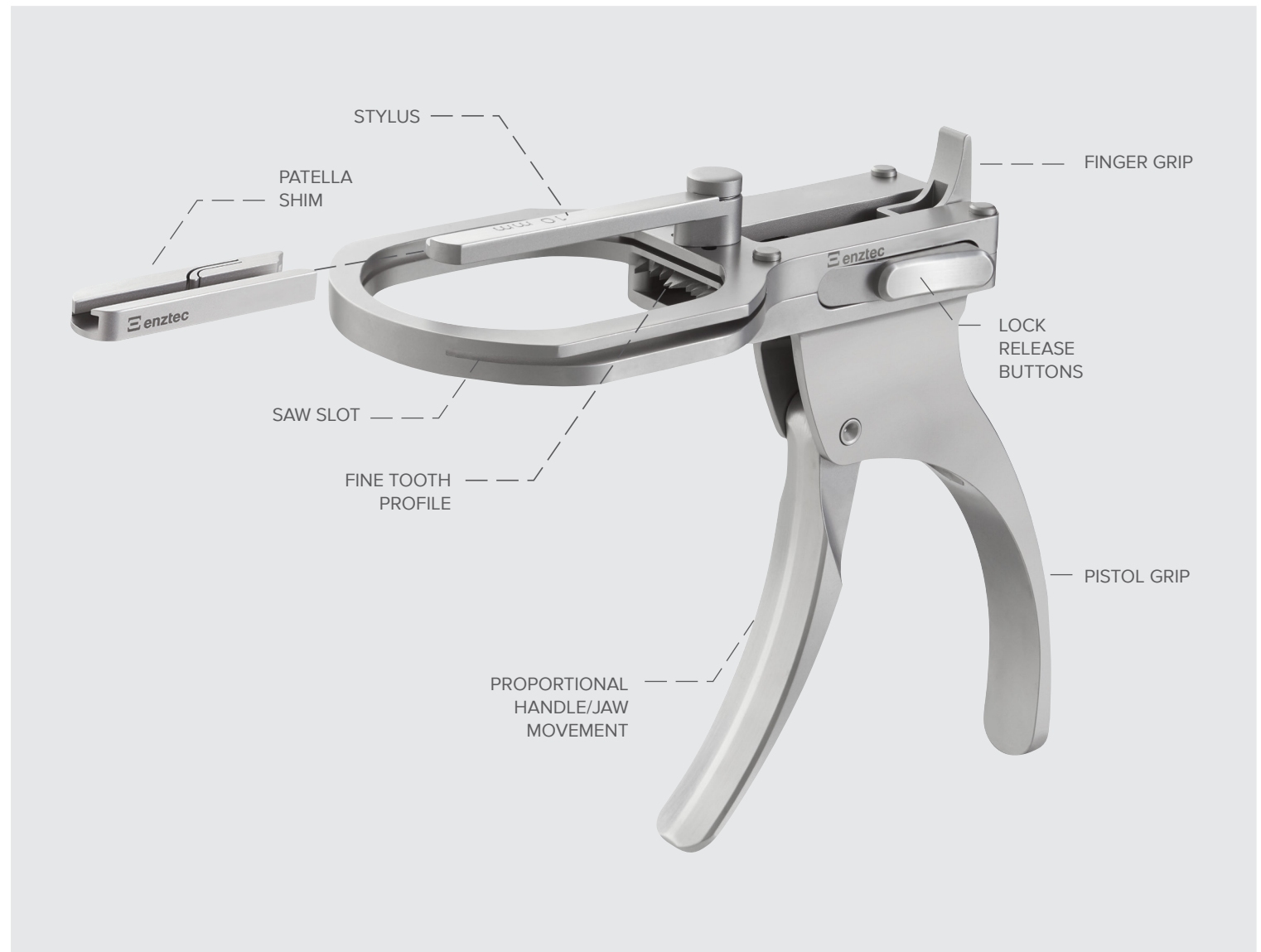

Enztec is focused on the design and manufacture of orthopaedic instruments. We are a proven, quality approved supplier to all major orthopaedic implant companies and welcome your enquiries.

For further information or to order your Enztec Premium Patella Saw Guide, visit: www.enztec.com
enquiries@enztec.com
NZ +64 3 348 0203

Enztec