

Classic Offset Cup Impactor+



# Classic Offset Cup Impactor+

The Enztec Classic OCI+ gives you the ability to check and adjust the position of the cup in-situ, giving you maximum placement control.

## Features:


- Detachable trinket allows for in-situ visualisation and intraoperative imaging of the cup before final impaction
- Infinite adjustment of cup position allows you to secure the cup in desired orientation prior to impaction
- Generous offset provides ample clearance of bone and soft tissues
- Suitable for lateral and direct anterior approaches
- Compatible with thread-interfaced cups, enquire about compatibility with your cup

## Specifications

Materials: Stainless Steel, Propylux HS  
 Component compatibility: Thread-interfaced cups  
 Weight: 1.3kg / 2.9lbs  
 Length: 420mm / 16.5in  
 Offset: 55mm / 2.2in  
 Add your logo: available



FDA  
COMPLIANT

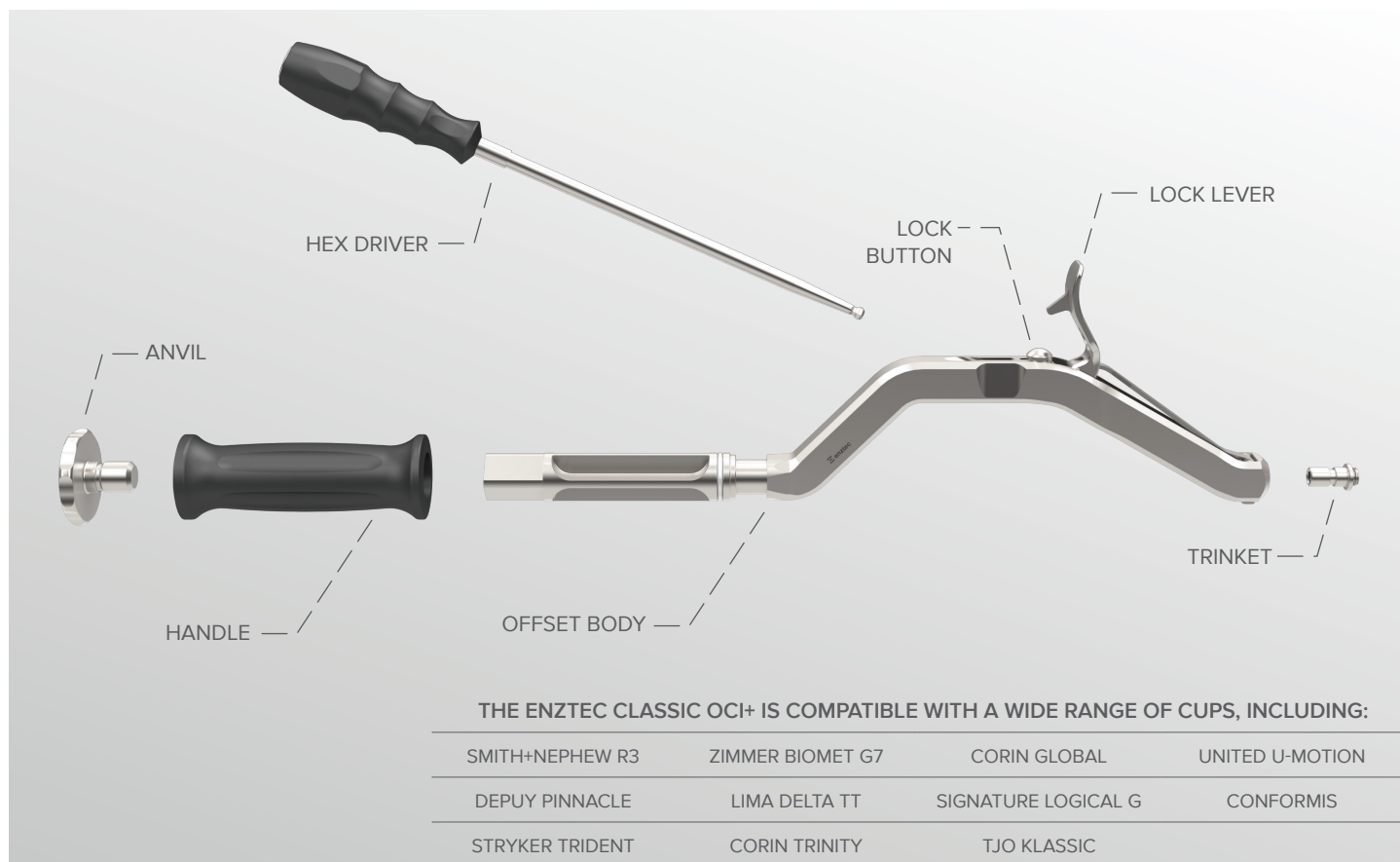

## THE ENZTEC CLASSIC OCI+ IS COMPATIBLE WITH A WIDE RANGE OF CUPS, INCLUDING:

SMITH+NEPHEW R3	ZIMMER BIOMET G7	CORIN GLOBAL	UNITED U-MOTION
DEPUY PINNACLE	LIMA DELTA TT	SIGNATURE LOGICAL G	CONFORMIS
STRYKER TRIDENT	CORIN TRINITY	TJO KLASSIC	

For further information or to order your Classic Offset Cup Impactor+, visit:

[www.enztec.com](http://www.enztec.com)

Enztec is an MDR-certified, ISO 13485 designer and manufacturer of orthopaedic instruments. We are a proven, quality approved supplier to all major orthopaedic implant companies and welcome your enquiries.

[enquiries@enztec.com](mailto:enquiries@enztec.com) | NZ +64 3 348 0203

